



Filter Frequency Response Characterization

Tool Reference

RST Reference Number: RST24CV14.01

Date of Publication: 02/07/2024

Recommended Citation: U.S. Food and Drug Administration. (2024). *Mock Circulatory Loop to Generate Variable Adult Heart Conditions for Evaluating Mechanical Circulatory Support Devices* (RST24CV14.01). <https://cdrh-rst.fda.gov/mock-circulatory-loop-generate-variable-adult-heart-conditions-evaluating-mechanical-circulatory>

For more information

[Catalog of Regulatory Science Tools to Help Assess New Medical Devices](#)

Disclaimer

About the Catalog of Regulatory Science Tools

The enclosed tool is part of the Catalog of Regulatory Science Tools, which provides a peer-reviewed resource for stakeholders to use where standards and qualified Medical Device Development Tools (MDDTs) do not yet exist. These tools do not replace FDA-recognized standards or MDDTs. This catalog collates a variety of regulatory science tools that the FDA's Center for Devices and Radiological Health's (CDRH) Office of Science and Engineering Labs (OSEL) developed. These tools use the most innovative science to support medical device development and patient access to safe and effective medical devices. If you are considering using a tool from this catalog in your marketing submissions, note that these tools have not been qualified as [Medical Device Development Tools](#) and the FDA has not evaluated the suitability of these tools within any specific context of use. You may [request feedback or meetings for medical device submissions](#) as part of the Q-Submission Program.


For more information about the Catalog of Regulatory Science Tools, email RST_CDRH@fda.hhs.gov.

Filter Frequency Response Characterization

***Note:** The test fluid used is blood analog fluid (BAF). Please prepare the BAF by following detailed instructions outlined in the document, '[Preparation of Blood Analog Fluid \(BAF\)](#)'. For preliminary or feasibility testing, De-ionized (DI) water may be used in place of BAF as the test fluid. While all test protocols in this document are written for BAF, they are also directly applicable to DI water, when used as the test fluid.*

1. Connect the simplified flow loop per the diagram shown below.

***Note:** The diagram below shows representative samples of 1 pressure transducer and 1 flow probe. The flow loop may be modified to incorporate 3 pressure transducers and 2 flow probes to achieve simultaneous sensor characterization, if needed.*

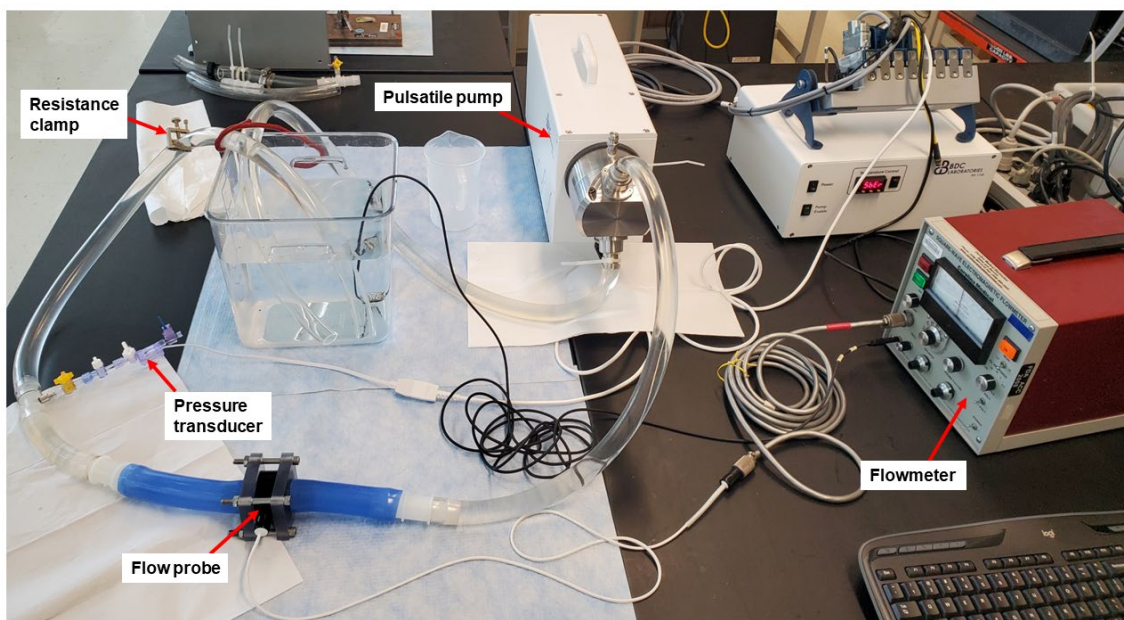

Available pressure transducers and amplifier:

Three (left ventricle (LV), left atrium (LA), and aorta (Ao)) pressure transducers [Utah Medical Products Inc., Midvale, Utah] connected to the pressure amplifier module [Bus 21097, ViVitro Labs, Inc., Victoria, BC, Canada].

Available flow sensors and flowmeter modules:

Flowmeter	Flow sensor
T402 (dual channel), Transonic Systems, Inc., Ithaca, NY	ME-25PXN, ME-11PXL, ME-16PXL
FM501, Carolina Medical Electronics, East Bend, NC	EP688 (size 88 mm cir)

***Note:** Since the EP688 flow probe requires a conductive fluid, sodium chloride (NaCl) would need to be mixed with DI water, when used as the test fluid.*


2. Make the appropriate electrical connections:
 - a. Connect the flow sensor/pressure transducer to the corresponding flowmeter/pressure amplifier module.
 - b. Connect the pulsatile pump [PD-1100, BDC Laboratories, Wheat Ridge, CO] per its user manual.
 - c. Connect the output channels on the flowmeter/pressure amplifier module to the channels on the analog input (AI) module [NI-9205, National Instruments, Austin, TX] of the data acquisition system [cDAQ-9174, National Instruments, Austin, TX] per the table below.

Instrument Channel	AI Channel on cDAQ-9174
Pressure transducer amplifier – Aorta (Ao) pressure channel	AI 1, NI-9205
Pressure transducer amplifier – Left ventricle (LV) pressure channel	AI 2, NI-9205
Pressure transducer amplifier – Left Atrium (LA) pressure channel	AI 3, NI-9205
T402 output channel	AI 4, NI-9205
FM501 output channel	AI 5, NI-9205

- d. Connect the cDAQ data cable to the mock circulatory loop (MCL) laptop.

***Note:** The pulsatile pump [PD-1100, BDC Laboratories, Wheat Ridge, CO] is controlled using its own desktop computer while the MCL laptop and data acquisition system is used to monitor and acquire data from pressure transducers and flow probes/sensors.*


***Note:** As an alternative to the MCL laptop and data acquisition, the pulsatile pump's [PD-1100, BDC Laboratories, Wheat Ridge, CO] data monitoring and acquisition capabilities can be used which allows for 1 pressure transducer and 1 flow probe at any instant of time. Please refer to the PD-1100 user manual.*
3. Power ON all equipment.
4. Open and Run the LabVIEW VI (or Virtual Instrument) that executes the custom MCL software (snapshot graphical user interface (GUI) shown below).

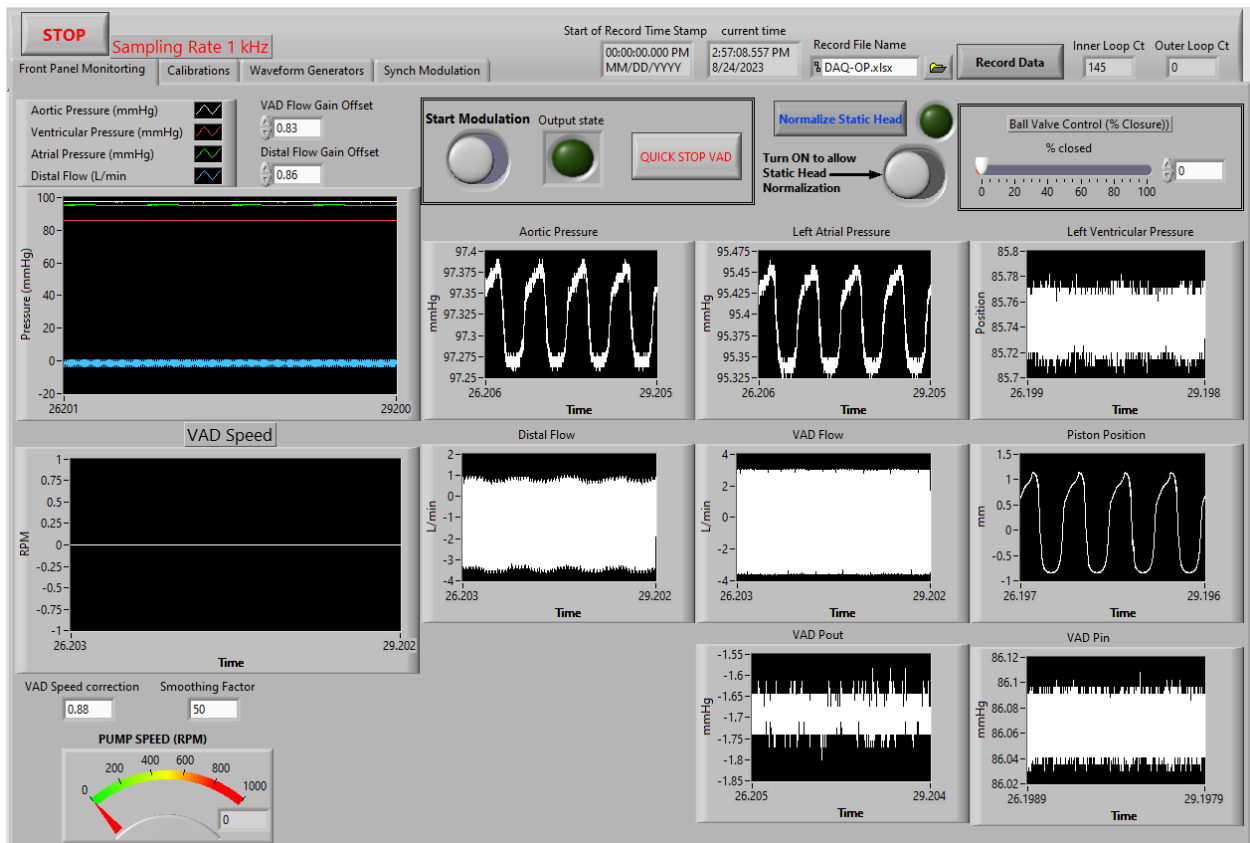


5. Apply the sensor/transducer Calibration Slope and Calibration Intercept such that the MCL software displays the pressure in mmHg and flow rate in L/min. Please refer to the documents, '[Pressure Transducer Calibration](#)' and '[Flow Sensor Calibration](#)'.


Note: The MCL software has been developed with the goal of evaluating ventricular assist devices (VADs) within the MCL, which is the reason for the nomenclatures, 'Distal flow' and 'VAD flow'. The MCL is designed to incorporate a maximum of 2 flow sensors at any given time of testing. Thus, in the absence of a VAD connected to the MCL, the 'VAD flow' AI channel is used to record cardiac output (CO) using the EP688 flow probe. 'Distal flow' represents the aortic (Ao) or systemic flow recorded by the ME-25PXN flow sensor.

Note: Please ensure that the VAD Flow Gain Offset and Distal Flow Gain Offset (in the Front Panel Monitoring tab) are set to 0.00.

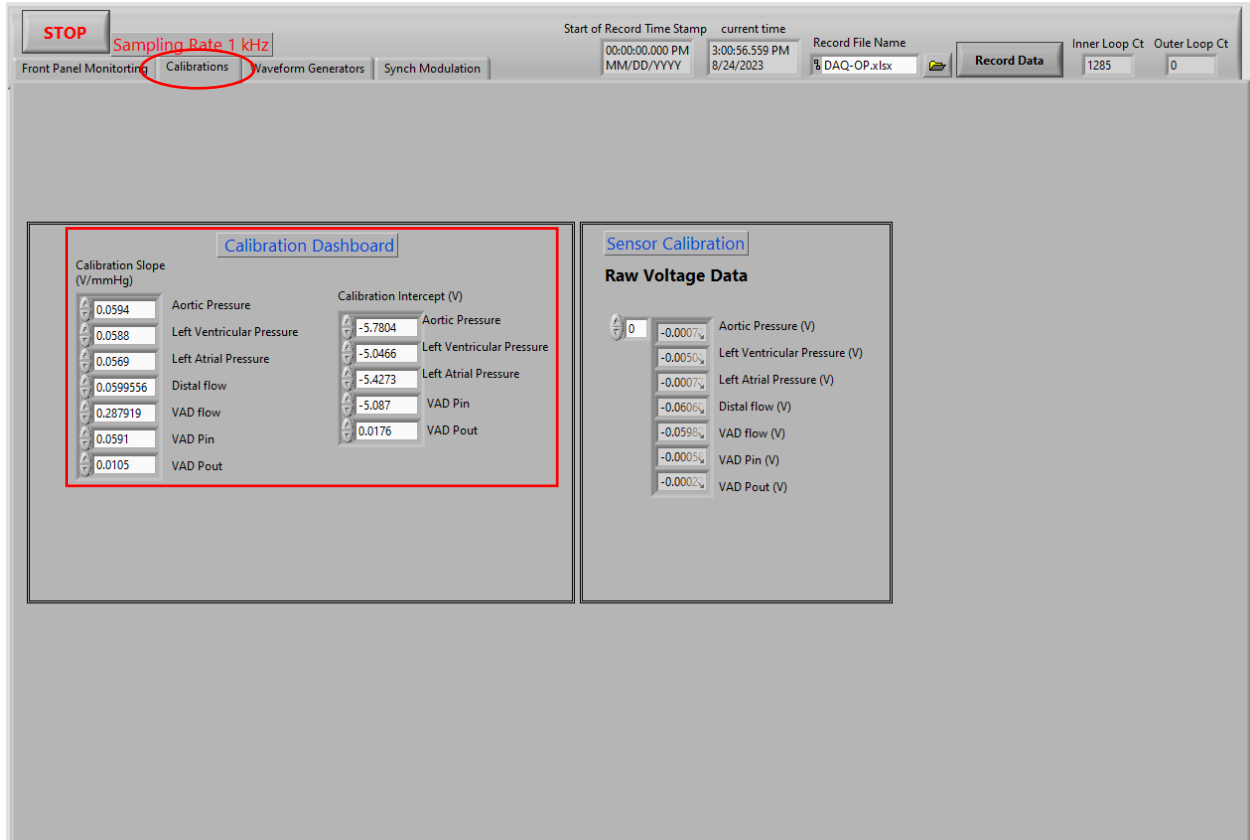

6. Run the pulsatile pump [PD-1100, BDC Laboratories, Wheat Ridge, CO] per its user manual using its accompanying software [Statys® PD, BDC Laboratories, Wheat Ridge, CO]. Use the following Flow Controller > Pump Control settings: Rate: 110 BPM, Target Flow: 3.5 LPM, Profile: 35% Systole.
7. Start to close the resistance valve till the peak pressure is under 300 mmHg and the flow rate is ~2.8 L/min – 3 L/min. A pressure pulse reading of ~ 0/250 mmHg is acceptable.
8. Set the lowest low pass filter /pulsatile frequency response setting on the flowmeter/pressure amplifier module.
9. Using the MCL software, record ~ 10 seconds of data by selecting Record Data (Note: A green illuminated Recording tab will appear when the software is recording data). Select an appropriate file name and path under Record file Name to save the Excel (.xlsx) data file. At the end of 10 seconds, select the green illuminated Recording tab to stop recording. Please refer to the snapshots below.

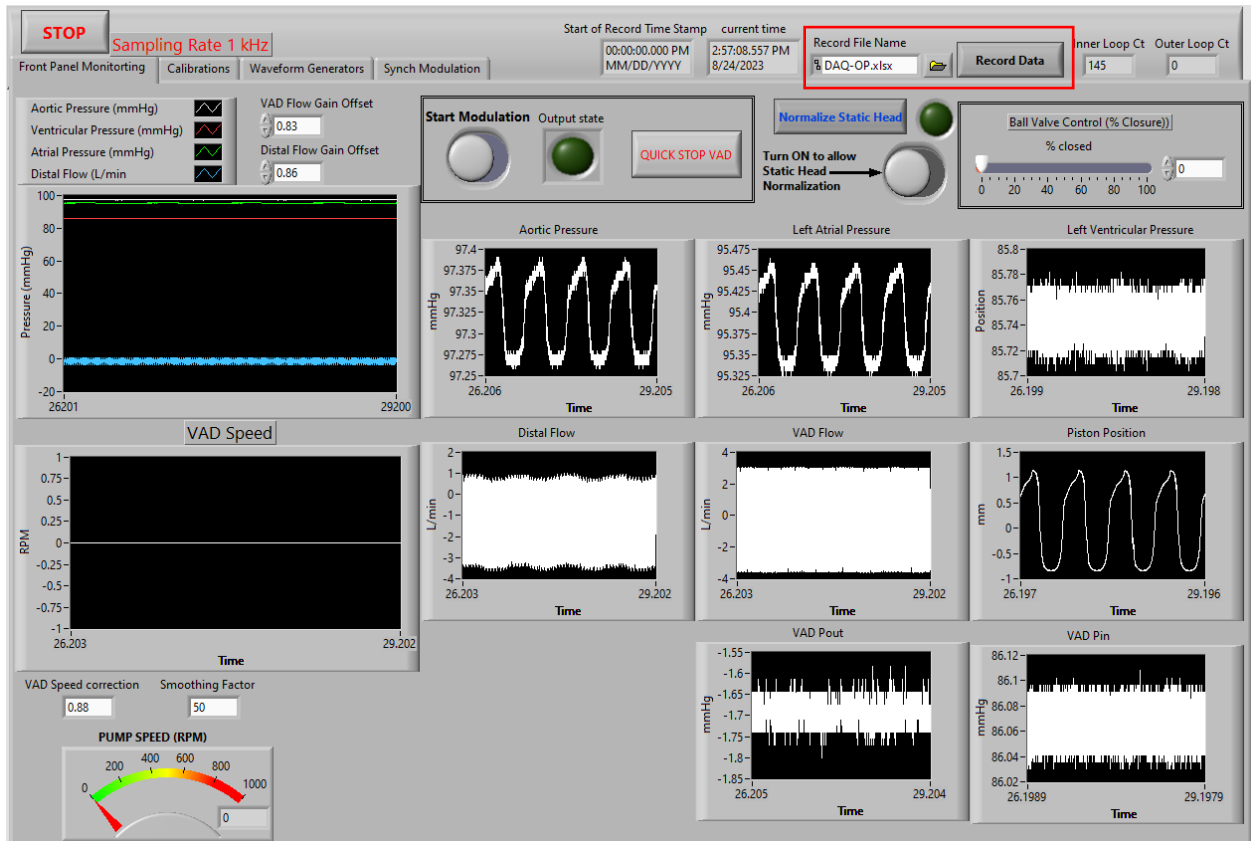

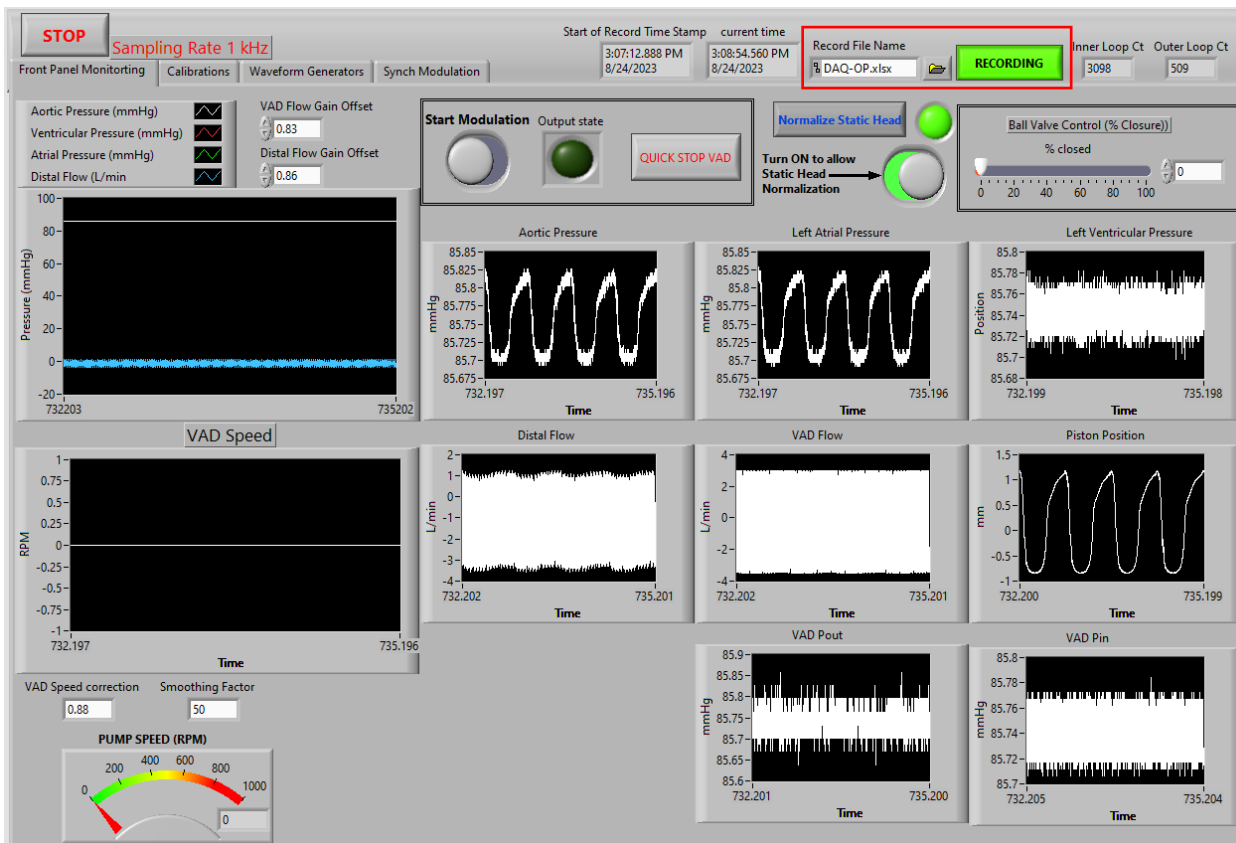

10. Repeat step 10 at all other low pass filter/pulsatile frequency response settings.
11. Plot flow rate/pressure vs. time over one beat cycle for all the low pass filter/pulsatile frequency response settings. Record the date of filter frequency response characterization.